

**"FINALLY A WHEELCHAIR THAT GETS ME TO THE TOP OF MY DRIVEWAY. SOMETIMES IT'S THE LITTLE THINGS LIKE BEING ABLE TO COLLECT MY OWN MAIL THAT MAKE A DIFFERENCE"**

**Sue**  
United Kingdom



**magic**  
*Mobility*<sup>®</sup>

Frontier V4 FWD  
**FOR THE ROAD LESS TRAVELLED**

### FULL LIGHTING PACKAGES

Been out having too much fun and it's getting late? Put your lights on.



### TAKE CONTROL

An extensive range of advanced control systems and environmental controls are available. Easy to program and use.



### REHAB SEATING

Infinite options for infinite body shapes. Our rehab seating range is easily configurable to suit your individual needs.

### TAILORED FOR YOU

Select the right headrest, backrest, armrests, legrests, seating and cushioning to suit your lifestyle.



CE UK  
CA

V4/AUST/AUGY24

Frontier  
**V4**<sup>®</sup>





# FRONTIER V4 FRONT WHEEL DRIVE

A powered wheelchair for people who know what they want. The V4 FWD offers a great solution to some of life's most tricky corners. The thinking person's chair, the V4 FWD is designed for those wanting to get around their homes with ease, using desks and tables, or driving frequently. Choose Off-road, Crossover or Urban wheels to ensure you get the best powerchair for your lifestyle.

"WE'VE BEEN  
CREATING  
HIGH-QUALITY  
WHEELCHAIRS  
FOR 30 YEARS,  
SO THAT YOU  
CAN LIVE THE  
LIFE YOU  
CHOOSE."

Ashley Daff  
R&D Director




## ATTENTION TO DETAIL

We know that small touches make a big difference. Our powerchair range is configurable and modifiable to suit your individual needs.

## ROBUST QUALITY

We never compromise on quality. Your powerchair is your mobility, therefore reliability is paramount.

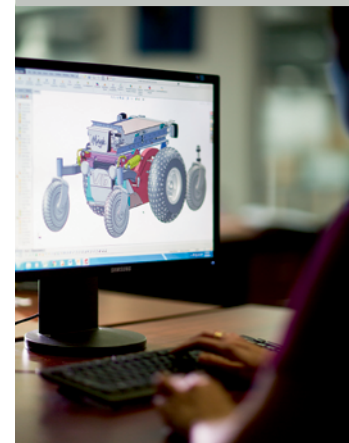

## ENGINEERED TO THE LAST BOLT

We listen to and value your feedback. Our team of in-house engineers work side by side with end users and therapists. We're all about innovation; if we think we can make our powerchairs better, we will.



## AUSTRALIAN DESIGNED AND MADE

We are proud that our powerchair range is designed and made in our Australian manufacturing facility.

### Basic specifications:\*

- Length: 995mm
- Width: 710mm - AT wheels  
700mm - Crossover wheels  
640mm - Urban wheels
- Seat to floor height: 435mm
- Base turning radius: 520mm
- Weight limit: 182kg (400lbs),  
155kg (340lbs) with lift and/or tilt
- Top speed: 10kph (6.2mph)
- Motors: 800W Gear in Line
- Wheels: Off-road 14" black,  
Crossover 14" black or grey,  
Urban 14"x3" black or grey  
(solid or pneumatic)
- Castors: 220mm (2.80/2.50-4),  
black or grey (solid or pneumatic)
- Gel batteries
- Penny and Giles controller with  
LCD display (R-Net)
- Centre post legrest with one piece  
footplate
- Rehab seat pan (from 12"-24" wide)
- Crash test compliant tie down  
points

### Rehab positioning:

- Power options: 300mm (12")  
seat lift, tilt (45 degrees), lift and  
tilt, transfer lift and tilt only, transfer  
tilt with lift and tilt, power recline  
and more
- Seating and backrests: MPS, rehab  
backrest, adjustable seat pans and  
seating
- Headrests and mounts: MPS,  
rehab
- Legrests: centre post, side mount,  
swingaway, power, gas or manual  
elevating
- Footplates: 1 or 2 piece fixed/flip  
up/angle adjustable, any specific  
requirement
- An extensive range of armpads,  
supports, restraints and pads are  
available

### The right wheelchair for:\*\*

- Adventurous off-road recreation
- Country living and working
- Urban outdoor living and working
- Roomy indoor spaces

### Options and accessories:

- One-click activator  
(CJSM 2 or JSM joystick)
- Fold forward back
- Retractable docking pin
- Puncture resistant tyre liners  
(Off-road and Crossover  
wheels only)
- 24V LED lighting kit
- 12V accessory socket and  
2 x USB charger
- Wheelchair jack
- Tray tables
- Luggage rack, bag hooks, bags,  
pockets
- Special controls: attendant,  
head, chin, foot, sip and puff,  
environmental controls
- RAM® mounting system (tablet  
and swing away tablet holder,  
mobile phone holder, swing away  
camera mount, attendant cup  
holder, action camera holder and  
cup holder)

### Base and tyre options:



SPACE BLUE BASE  
BLACK OFF-ROAD TYRES & RIMS



LUCENT LIME BASE  
BLACK CROSSOVER TYRES



BURNT ORANGE BASE  
BLACK URBAN TYRES

\*Specifications subject to change. Check online for current version  
\*\*Please remember that your wheelchair has electric motors and must  
not be driven through water, into rivers, creeks or the sea



Magic Mobility +61 3 8791 5600  
enquiries@magicmobility.com.au  
Explore with us [magicmobility.com.au](http://magicmobility.com.au)



Scan for agent locations