



APEX
MOBILITY

KLAXON BROCHURE

apexmobility.com.au | 1300 212 192

The APEX Mobility Product Range

The APEX Mobility team continually research new products and global innovation, which is why we can provide the very best assistive technology equipment available in the Australian market today. Before any products are added to the APEX Mobility range, they are tested and trialed to make sure they meet all functional and quality requirements, plus regulatory standards for safety.

Products featured in this brochure are part of our dynamic range of quality assistive technology equipment. Most of these products are scripted/prescribed, which means they are customised to suit the specific needs of individual clients.



*Klaxon Mobility
- the art of moving*

Twist	4
Klick Power	6
Klick Race	7
Klick Monster	8
Tetra Handlebar Options	9
Linking System Klaxon® Main Accessories	10

Handbike models for everyday life, trips and adventures.

Wheelchair users can now experience a new level of freedom without compromising the manoeuvrability of their manual wheelchair. Klaxon Klick or Twist devices are small, light, compact, and suitable for travel and everyday use. The Innovative Linking System Klaxon® means it can connect to any hand bike in seconds, without the need for personal assistance, and the universal hub adapts to any fixed or folding wheelchair frames.

The Electronic Brake System is a unique feature. One touch of the brake lever is all it takes to stop the Klick or Twist, without having to use the mechanical brake. During the EBS braking, the energy produced will be recovered and used to recharge the battery. The mechanical brake can still be activated by applying a greater force on the brake lever. The two braking systems will then merge in a single braking action.





Twist



Specifications:

Material	Aluminium
Motor	36 V x 750 W brushless
Wheel	10" smooth tread tyre
Brakes	with Controller: Electronic Brake System with Handlebar: 150 mm Mechanical Brake
Driving speed	with Controller: Front drive: up to 6 km/h Rear drive: up to 10 km/h with Handlebar: up to 10 km/h
Battery Information	Lithium 36 V x 4 Ah (weight 1.5 kg) Charging time: approx. 2h
Range	approx. 20 km (at moderate speed, on flat terrain and with user weighing 70 kg)
Weight	7.5 kg (without battery)
Connection	Available options: Central Linking System Klaxon® Lateral Linking System Klaxon® Rear Linking System Klaxon®
Dimensions	Twist: W 15 cm x H 34 cm x D 31 cm with Handlebar: W 45 cm x H 94 cm x D 42 cm
Technologies	Downhill Speed Control (in Controller Mode), Intelligent Cruise Control (in Controller Mode), Electronic Brake System
Included accessories	Battery, Battery Charger, Frontal Light, Rear Light
Included services	Klaxon® App
Optionals	Controller, Handlebar, Additional battery

Discover the Twist

A game-changer. The Klaxon Twist is an innovative power-assist device for manual wheelchairs.

- Multi-configuration allowing you to connect via the front or rear
- Lightweight and manoeuvrable to accompany you on every journey
- Compact and elegant design



Several possibilities of use

Twist engaged on the front for a relaxed and safe outdoor ride. Choose your speed with the Controller and move with your hands free. Connect the handlebar and go fast in total comfort! The rear-mounted Twist is ideal for indoor use.

Quick and easy to connect

Twist can be connected and disconnected in seconds and in total autonomy. Twist is the solution for your independent mobility. Discover a new feeling of freedom without giving up the comfort of moving with your wheelchair.



Stylish, technological and wireless

Connect your driving system via Bluetooth. Control and select the pre-settings directly from the Klaxon® App. Nothing stands between you and your freedom! Twist is cool and elegant, perfect for moving around in style.

Safe and secure

Twist is made of high-quality materials and meets the requirements of innovation and technology. It is designed to accompany you at all times in complete safety and durability. Integrated LED lights to ensure visibility at night!



Lightweight and manoeuvrable

Twist is agile. It rotates 360° and is ideal for narrow spaces. The integrated ergonomic handle makes it easy to hold and transport. We take care of every detail.

Compact and transportable

Explore new destinations with Twist. Discover the pleasure of travelling like never before! Twist is compact and lightweight and can be carried and transported by car, train or plane.



Suitable for all types of terrain

Twist overcomes everyday obstacles with ease. Manages climbs efficiently. When installed in the front it lifts the wheelchair's front wheels and allows you to move on all types of terrain.

Linking System Klaxon®

The innovative linking system can be mounted on any type of wheelchair, rigid or foldable. Expand the functionality without sacrificing manoeuvrability and lightness.



Klick Power



Specifications:

Material	Aluminium
Motor	48 V x 750 W brushless, all wheel
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 Speed settings, walking pace, USB charging port
Dedicated Buttons	Adaptive Cruise Control, Electronic Brake System Reverse gear, Double mapping 'Sport' and 'Comfort'
Battery	Lithium 48 V x 5.8 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle
Brakes	Dual Lever Brake 160 mm with EBS - Electronic Brake System
Size	W 50 cm x H 90 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	8.5 kg (battery excluded)
Range	About 50 km (speed setting 1, ground level)
Included accessories	Battery 11 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard



Klick Race

Specifications:



Material	Aluminium
Motor	48 V x 1500 W brushless, all wheel
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 Speed settings, walking pace
Battery	Lithium 48 V x 11 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle
Brakes	Dual Lever Brake 160 mm with EBS - Electronic Brake System
Size	W 50 cm x H 95 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	11.4 kg (battery excluded)
Range	About 50 km (speed setting 1, level ground)
Included accessories	Battery 11 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard



Klick Monster

Specifications:



Material	Aluminium
Motor	48 V x 1500 W brushless
Wheel	20" tyre XL
Display	Multifunction LCD, 5 Speed settings, cruise control, walking pace
Battery	Lithium 48 V x 11 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle
Brakes	Dual Lever Brake 160 mm with EBS - Electronic Brake System
Size	W 50 cm x H 90 cm x D 50 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	15.0 kg (battery excluded)
Range	About 40 km (speed setting 1, level ground)
Included accessories	Battery 11 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard

Tetra Handlebar Options

Klick Power Tetra



Specifications:

Material	Aluminium
Motor	48 V x 750 W brushless
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 Speed settings, walking pace, USB charging port
Battery	Lithium 48 V x 5.8 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle; Comfort Drive System (accelerator and brake activated through push and pull movement)
Brakes	Double Brake System 160 mm
Size	W 50 cm x H 95 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	9.1 kg (battery excluded)
Range	About 50 km (speed setting 1, ground level)
Included accessories	LCD display, bell, foldable kickstand, splashguard

Klick Race Tetra



Specifications:

Material	Aluminium
Motor	48 V x 1500 W brushless, all wheel
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 Speed settings, walking pace
Battery	Lithium 48 V x 11 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle; Comfort Drive System (accelerator and brake activated through push and pull movement)
Brakes	Double Brake System 160 mm
Size	W 50 cm x H 95 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	12.0 kg (battery excluded)
Range	About 50 km (speed setting 1, ground level)
Included accessories	Comfort Drive System, Battery 11 Ah with battery charger, LCD display, bell, foldable kickstand, splashguard

Linking System Klaxon®

Worldwide Patented

The innovative linking system connects to any hand bike of the KLICK product range in seconds, allowing you to switch between your favourite Klick hand bikes with ease. The universal connection hub adapts to both folding and rigid wheelchair frames. Experience a new level of freedom without compromising the manoeuvrability and comfort of your manual wheelchair. LINE-UP. ENGAGE. KLICK.



QUICK AND EASY TO ATTACH



LINKS IN SECONDS



UNIVERSAL CONNECTION HUB

Tetra Linking Bar Option

Main Accessories



Mini Flybag



Standard Flybag



Race Battery (48V x 11Ah)
Fly-Version (48V x 5.8Ah)



Standard Light



Front Light Monster/Race




Additional Traction Weight



Phone Holder



X-Large Wheels Set

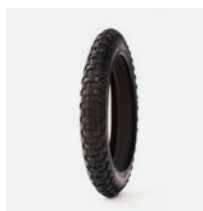

14" Tyres



Race Battery Charger




Power Battery Charger
Hybrid Battery Charger



Slim Battery Charger
Mini Battery Charger



Handlebar Extension Std
Handlebar Extension Tetra

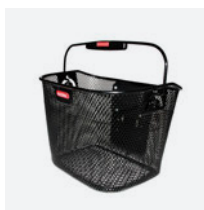


Basket
10L and 16L



Basket Linking Kit



Rear-view Mirror

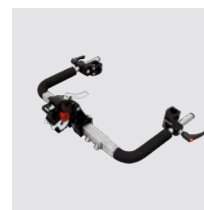

Lateral Linking System
Klaxon®



Central Linking System
Klaxon®



Product Demonstrations & Trials

APEX Mobility does not have demo units available. In many cases, Dealers are offered subsidised pricing for units that they will hold for showroom display or trial use only.

Dealers of Assistive Technology

Contact your APEX Mobility account manager to discuss your needs.

Therapists

To arrange a trial of any of this range, contact the following dealer:

Authorised Dealer



 1300 212 192  sales@apexmobility.com.au

 apexmobility.com.au

V2 May 2023